




## HAMMERTOE TRUSS SYSTEM™

### TITANIUM PROXIMAL INTERPHALANGEAL (PIP) JOINT FUSION DEVICE

The Hammertoe Truss System™ (HTS) has an Advanced Structural Design that incorporates 4WEB Medical's proprietary Truss Implant Technology™. Under normal loading conditions the struts in the truss implant transfer strain to adjacent cellular material which stimulates a mechanobiologic response.

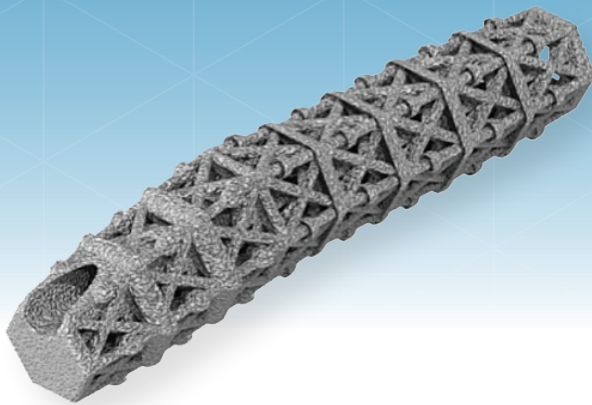




- ▲ Distribution of load through the implant struts delivers strain to adjacent cellular material which stimulates a mechanobiologic response.
- ▲ Open architecture design allows for greater graft volume and bone growth throughout the entire construct.\*
- ▲ Truss Implant design provides maximum strength with a minimal amount of material, which limits imaging artifacts.

Sterile-packed cannulated  
implant & drill bit for  
hospital efficiency and  
patient safety


Proximal and distal  
implant taper allows for a  
compression-fit

Multiple diameters and 0°  
and 10° options to  
accommodate varying  
patient anatomy

Single-use, sterile packed  
instrument kit designed  
for improved surgical  
efficiency

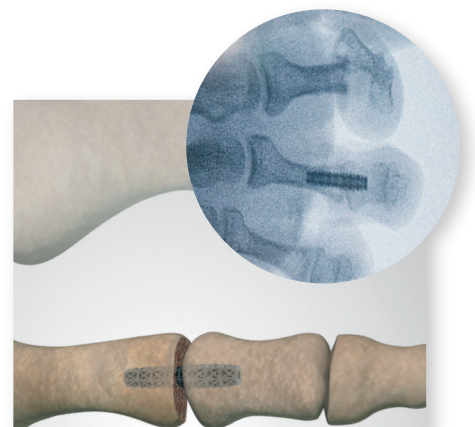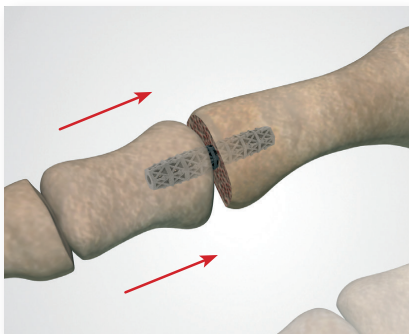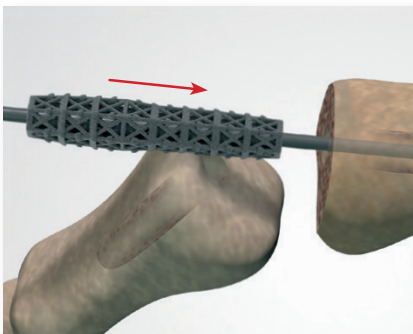

\*Data on File

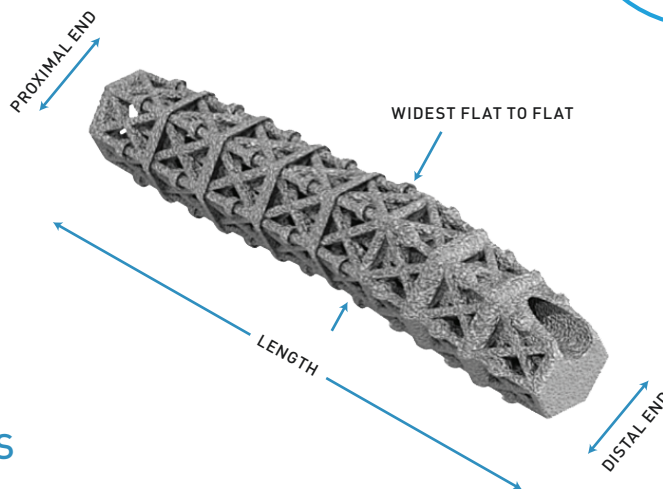

## HTS IMPLANT SPECIFICATIONS

### PRODUCT CATALOG

CATALOG #	DESCRIPTION	WIDEST FLAT TO FLAT	PROXIMAL END	DISTAL END	LENGTH	DEGREE
0 DEGREE*						
HTS-XS1500-SP	Hammertoe 3.0 x 15mm, 0°	3.0mm	2.3mm	2.4mm	15mm	0°
HTS-SM1700-SP	Hammertoe 3.2 x 17mm, 0°	3.2mm	2.5mm	2.6mm	17mm	0°
HTS-MD1900-SP	Hammertoe 3.5 x 19mm, 0°	3.5mm	2.6mm	2.8mm	19mm	0°
HTS-LG2100-SP	Hammertoe 3.8 x 21mm, 0°	3.8mm	2.8mm	3.0mm	21mm	0°
HTS-XL2600-SP	Hammertoe 4.7 x 26mm, 0°	4.7mm	3.7mm	4.0mm	26mm	0°
10 DEGREE*						
HTS-XS1510-SP	Hammertoe 3.0 x 15mm, 10°	3.0mm	2.3mm	2.4mm	15mm	10°
HTS-SM1710-SP	Hammertoe 3.2 x 17mm, 10°	3.2mm	2.5mm	2.6mm	17mm	10°
HTS-MD1910-SP	Hammertoe 3.5 x 19mm, 10°	3.5mm	2.6mm	2.8mm	19mm	10°
HTS-LG2110-SP	Hammertoe 3.8 x 21mm, 10°	3.8mm	2.8mm	3.0mm	21mm	10°
HTS-XL2610-SP	Hammertoe 4.7 x 26mm, 10°	4.7mm	3.7mm	4.0mm	26mm	10°

\*Each implant is packaged with the corresponding drill bit

## HTS INSTRUMENT KIT

HTS-KIT001-SP

CATALOG #	DESCRIPTION	QTY
HTS-000001	HTS Insertion Handle	1
HTS-000002	HTS Sizer	1
HTS-000003	HTS Tamp	1
HTS-000004	K-wire	3

