




LATERAL SPINE *TRUSS* SYSTEM™

4WEB's *TRUSS* Implant Technology™ redefined the standard of care as the first FDA-cleared 3D printed interbody fusion device with an Advanced Structural Design. The LSTS continues the legacy of delivering innovation, strength, and performance to every procedure.

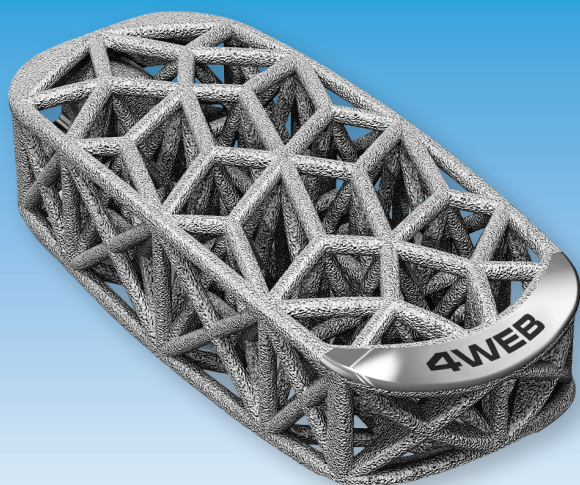

THERAPEUTIC BY DESIGN™




- ▶ Truss Implant Technology™ provides a synergistic relationship between the hierarchical surface roughness and the macro architecture.¹
- ▶ After surface features initiate osteogenic activity, the structural mechanics of the truss implant provide a mode of action that converts physiologic load into therapeutic strain at the surgical site, amplifying osteogenic activity in local osteo-biologic material.²
- ▶ This therapeutic effect provides ongoing support while biologic fusion takes place.³

Open Architecture Design allows for bone growth throughout the entire implant construct³

Snowshoe Interface distributes load across the bone contact surface, reducing the risk of subsidence³

Lattus360 Access System for use in both direct lateral and ATP procedures*

Bi-convex anatomic structure is designed to increase the implant-to-endplate contact surface

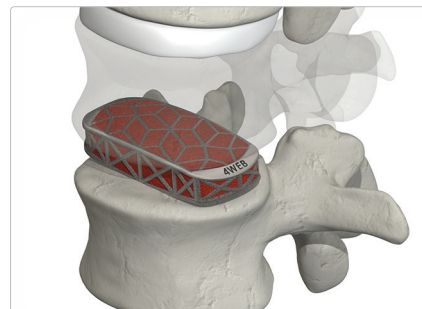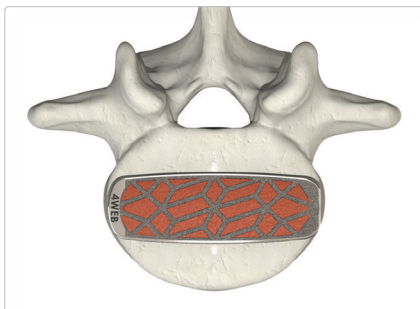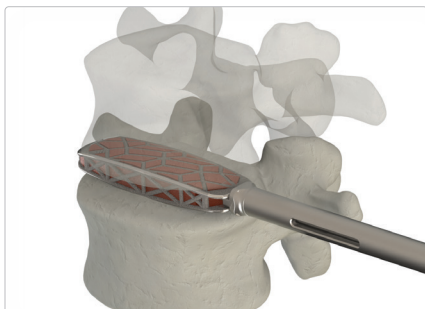

¹ Rowe, et al., SMISS, Annual Forum '19, p.52


² Lee et al. Adv. Healthcare Mater. 2025, 2501856

³ Data on file

* Lattus360 is commercially not available in the EU



VISIT
PRODUCT
PAGE

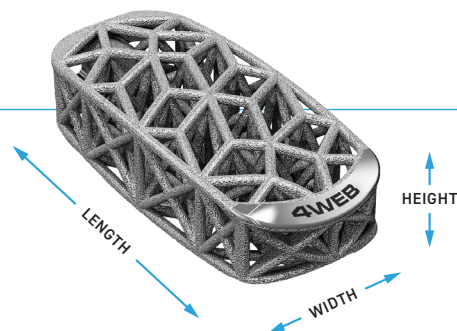

LSTS IMPLANT SPECIFICATIONS

CATALOG #	FOOTPRINT (W, L, H)	LORDOSIS
0 DEGREE		
LSTS-SM0008-45-SP	18 x 45 x 8mm	0°
LSTS-SM0010-45-SP	18 x 45 x 10mm	0°
LSTS-SM0012-45-SP	18 x 45 x 12mm	0°
LSTS-SM0014-45-SP	18 x 45 x 14mm	0°
LSTS-SM0008-50-SP	18 x 50 x 8mm	0°
LSTS-SM0010-50-SP	18 x 50 x 10mm	0°
LSTS-SM0012-50-SP	18 x 50 x 12mm	0°
LSTS-SM0014-50-SP	18 x 50 x 14mm	0°
LSTS-SM0008-55-SP	18 x 55 x 8mm	0°
LSTS-SM0010-55-SP	18 x 55 x 10mm	0°
LSTS-SM0012-55-SP	18 x 55 x 12mm	0°
LSTS-SM0014-55-SP	18 x 55 x 14mm	0°
LSTS-MD0008-45-SP	22 x 45 x 8mm	0°
LSTS-MD0010-45-SP	22 x 45 x 10mm	0°
LSTS-MD0012-45-SP	22 x 45 x 12mm	0°
LSTS-MD0014-45-SP	22 x 45 x 14mm	0°
LSTS-MD0008-50-SP	22 x 50 x 8mm	0°
LSTS-MD0010-50-SP	22 x 50 x 10mm	0°
LSTS-MD0012-50-SP	22 x 50 x 12mm	0°
LSTS-MD0014-50-SP	22 x 50 x 14mm	0°
LSTS-MD0008-55-SP	22 x 55 x 8mm	0°
LSTS-MD0010-55-SP	22 x 55 x 10mm	0°
LSTS-MD0012-55-SP	22 x 55 x 12mm	0°
LSTS-MD0014-55-SP	22 x 55 x 14mm	0°
6 DEGREE		
LSTS-SM0608-45-SP	18 x 45 x 8mm	6°
LSTS-SM0610-45-SP	18 x 45 x 10mm	6°
LSTS-SM0612-45-SP	18 x 45 x 12mm	6°
LSTS-SM0614-45-SP	18 x 45 x 14mm	6°
LSTS-SM0608-50-SP	18 x 50 x 8mm	6°
LSTS-SM0610-50-SP	18 x 50 x 10mm	6°
LSTS-SM0612-50-SP	18 x 50 x 12mm	6°
LSTS-SM0614-50-SP	18 x 50 x 14mm	6°
LSTS-SM0608-55-SP	18 x 55 x 8mm	6°
LSTS-SM0610-55-SP	18 x 55 x 10mm	6°
LSTS-SM0612-55-SP	18 x 55 x 12mm	6°
LSTS-SM0614-55-SP	18 x 55 x 14mm	6°
LSTS-SM0608-60-SP	18 x 60 x 8mm	6°
LSTS-SM0610-60-SP	18 x 60 x 10mm	6°
LSTS-SM0612-60-SP	18 x 60 x 12mm	6°
LSTS-SM0614-60-SP	18 x 60 x 14mm	6°
LSTS-MD0608-45-SP	22 x 45 x 8mm	6°
LSTS-MD0610-45-SP	22 x 45 x 10mm	6°
LSTS-MD0612-45-SP	22 x 45 x 12mm	6°
LSTS-MD0614-45-SP	22 x 45 x 14mm	6°

CATALOG #	FOOTPRINT (W, L, H)	LORDOSIS
LSTS-MD0608-50-SP	22 x 50 x 8mm	6°
LSTS-MD0610-50-SP	22 x 50 x 10mm	6°
LSTS-MD0612-50-SP	22 x 50 x 12mm	6°
LSTS-MD0614-50-SP	22 x 50 x 14mm	6°
LSTS-MD0608-55-SP	22 x 55 x 8mm	6°
LSTS-MD0610-55-SP	22 x 55 x 10mm	6°
LSTS-MD0612-55-SP	22 x 55 x 12mm	6°
LSTS-MD0614-55-SP	22 x 55 x 14mm	6°
LSTS-MD0608-60-SP	22 x 60 x 8mm	6°
LSTS-MD0610-60-SP	22 x 60 x 10mm	6°
LSTS-MD0612-60-SP	22 x 60 x 12mm	6°
LSTS-MD0614-60-SP	22 x 60 x 14mm	6°
12 DEGREE		
LSTS-SM1208-45-SP	18 x 45 x 8mm	12°
LSTS-SM1210-45-SP	18 x 45 x 10mm	12°
LSTS-SM1212-45-SP	18 x 45 x 12mm	12°
LSTS-SM1214-45-SP	18 x 45 x 14mm	12°
LSTS-SM1208-50-SP	18 x 50 x 8mm	12°
LSTS-SM1210-50-SP	18 x 50 x 10mm	12°
LSTS-SM1212-50-SP	18 x 50 x 12mm	12°
LSTS-SM1214-50-SP	18 x 50 x 14mm	12°
LSTS-SM1208-55-SP	18 x 55 x 8mm	12°
LSTS-SM1210-55-SP	18 x 55 x 10mm	12°
LSTS-SM1212-55-SP	18 x 55 x 12mm	12°
LSTS-SM1214-55-SP	18 x 55 x 14mm	12°
LSTS-SM1208-60-SP	18 x 60 x 8mm	12°
LSTS-SM1210-60-SP	18 x 60 x 10mm	12°
LSTS-SM1212-60-SP	18 x 60 x 12mm	12°
LSTS-SM1214-60-SP	18 x 60 x 14mm	12°
LSTS-MD1208-45-SP	22 x 45 x 8mm	12°
LSTS-MD1210-45-SP	22 x 45 x 10mm	12°
LSTS-MD1212-45-SP	22 x 45 x 12mm	12°
LSTS-MD1214-45-SP	22 x 45 x 14mm	12°
LSTS-MD1208-50-SP	22 x 50 x 8mm	12°
LSTS-MD1210-50-SP	22 x 50 x 10mm	12°
LSTS-MD1212-50-SP	22 x 50 x 12mm	12°
LSTS-MD1214-50-SP	22 x 50 x 14mm	12°
LSTS-MD1208-55-SP	22 x 55 x 8mm	12°
LSTS-MD1210-55-SP	22 x 55 x 10mm	12°
LSTS-MD1212-55-SP	22 x 55 x 12mm	12°
LSTS-MD1214-55-SP	22 x 55 x 14mm	12°
LSTS-MD1208-60-SP	22 x 60 x 8mm	12°
LSTS-MD1210-60-SP	22 x 60 x 10mm	12°
LSTS-MD1212-60-SP	22 x 60 x 12mm	12°
LSTS-MD1214-60-SP	22 x 60 x 14mm	12°

Note: Anterior face of implants are laser-marked with letter "A"

Hyperlordotic Implants available upon request